

Operator's Manual – Demo Mode

Dräger Infinity M540 Patient Monitor

w/Dräger Infinity C700 Monitor (Mirror)




Instructions for turning on the Dräger Infinity M540 patient monitor and Dräger Infinity C700 monitor (mirror), setting the Dräger Infinity M540 patient monitor to DEMO mode and silencing the ALARM warning. Also contains instructions for turning off the Infinity M540 and Infinity C700 monitors.

Dräger Infinity M540 Operation Manual for Demo Mode

Step 1 – Push In and Hold Green Power Button on Infinity M540 Until device begins to boot.



Step 2 – Push In and Hold Green Power Button on C700 Monitor until device begins to boot.

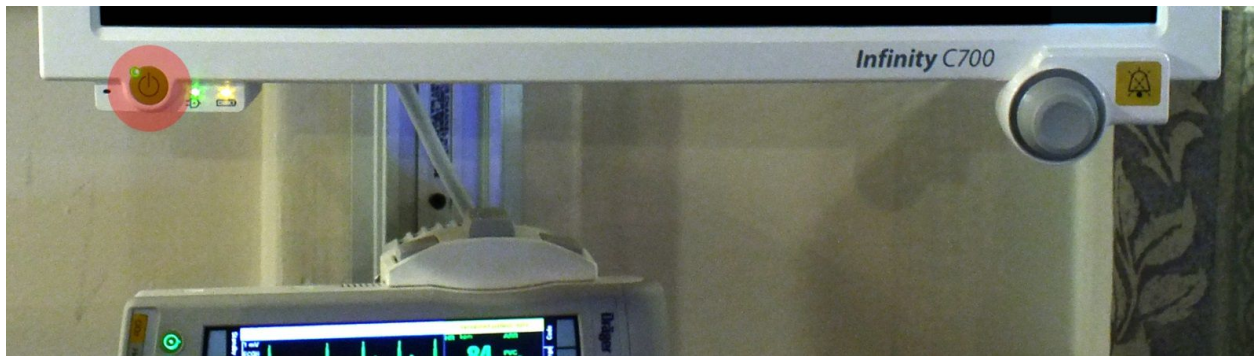


Step 3 – Press “Menu” Tab on touchscreen.




Dräger Infinity M540 Operation Manual for Demo Mode


Step 4 – Press “System Setup” tab on touchscreen.

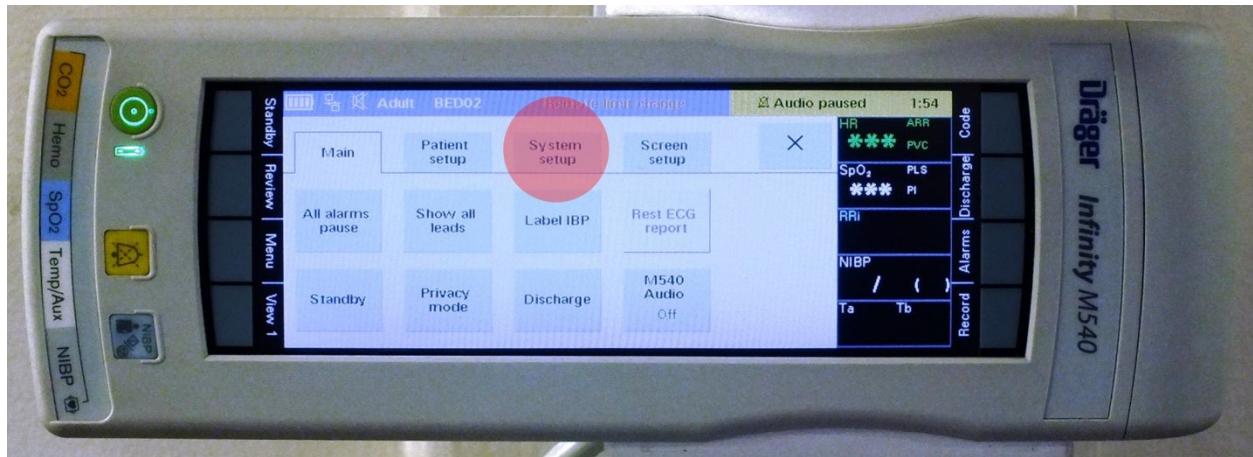

Step 5 – Press “Biomed” button on touchscreen.

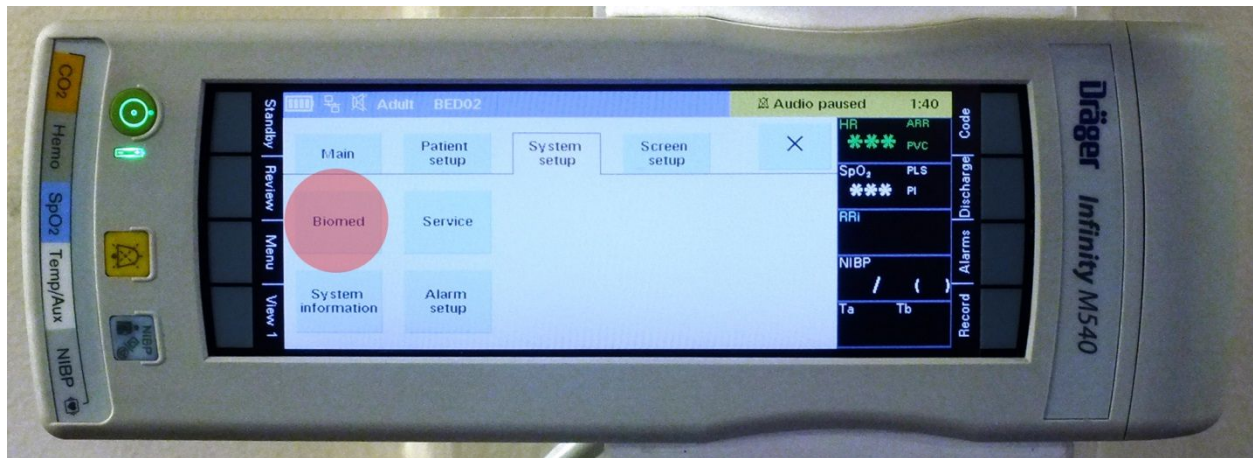


Step 6 – Press in sequence “4 – 7 – 5 – 6” on touchscreen.

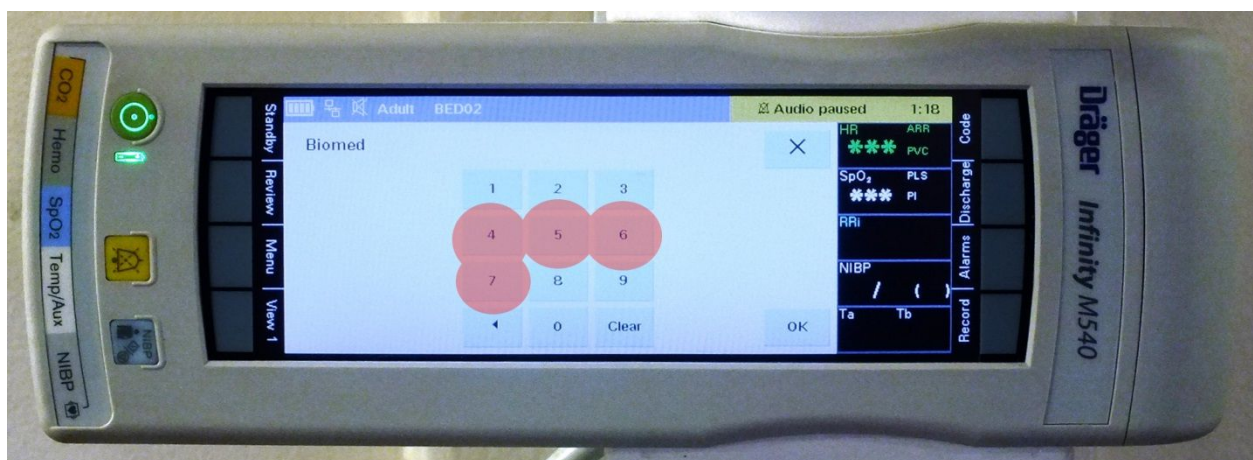


Dräger Infinity M540 Operation Manual for Demo Mode


Step 7 – Press “OK” on touchscreen.

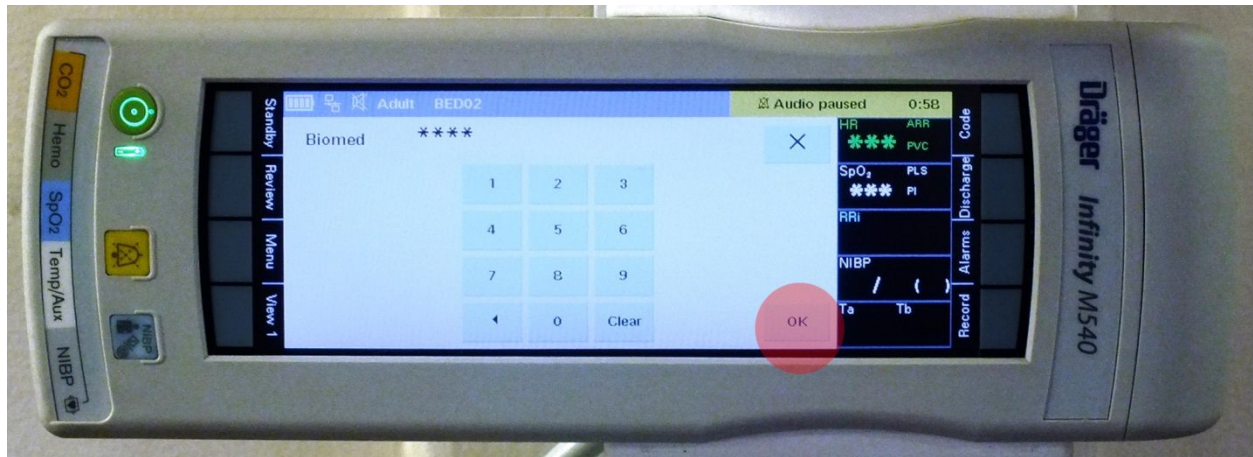

Step 8 – Press “Simulation” button on touchscreen.

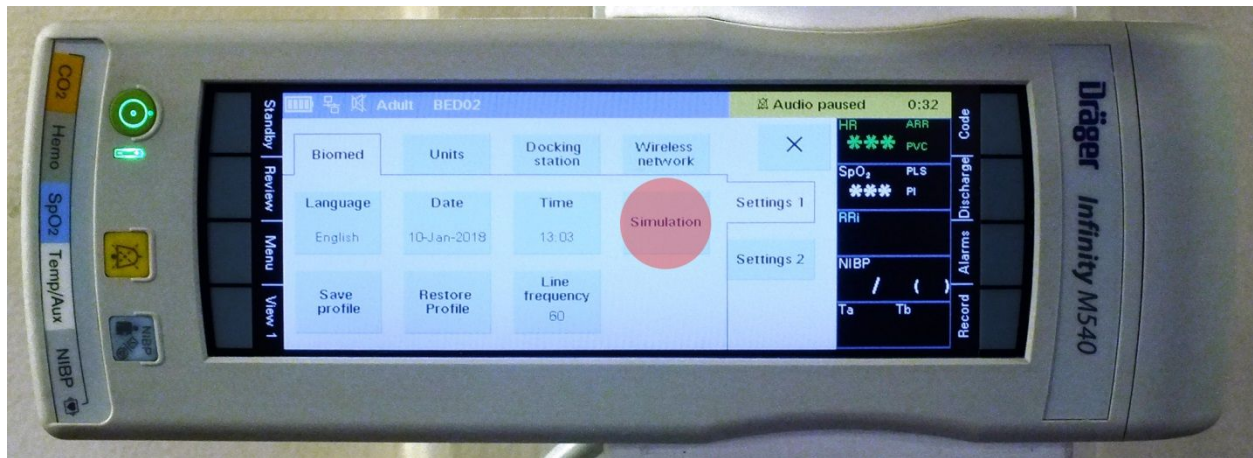


Step 9 – Press “OK” button on touchscreen.

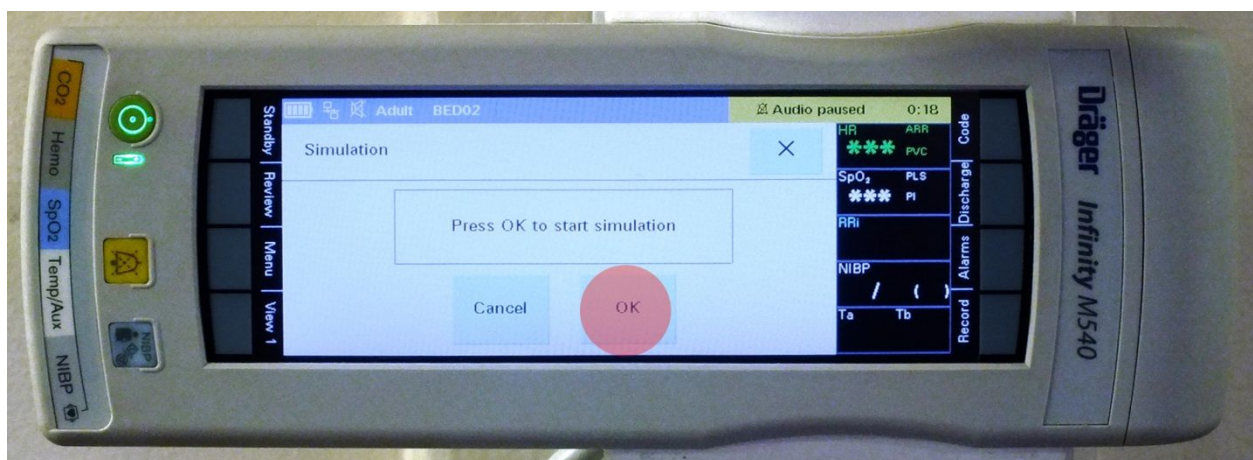

Dräger Infinity M540 Operation Manual for Demo Mode


Step 10 – Press Yellow Alarm button on Infinity C700 Monitor to temporarily silence alarm.

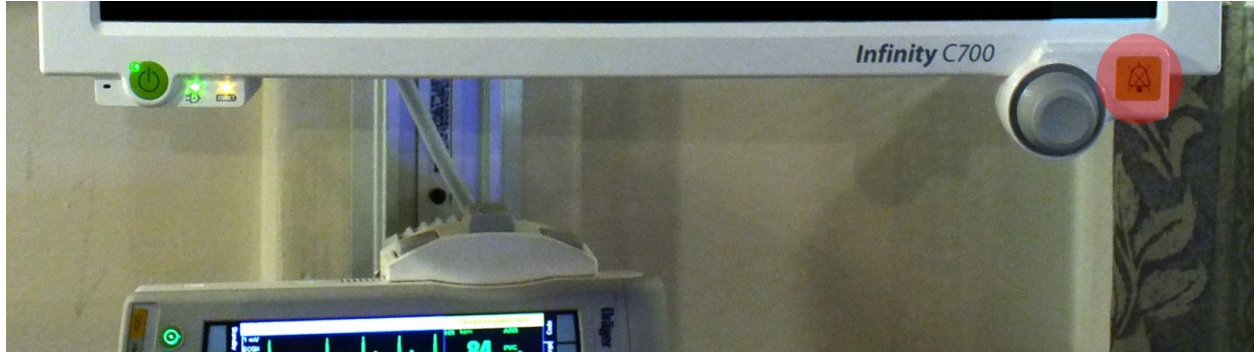

Step 11 – When alarm sounds press “84” on touchscreen.




Step 12 – Press upper range temp “75” button on touchscreen.

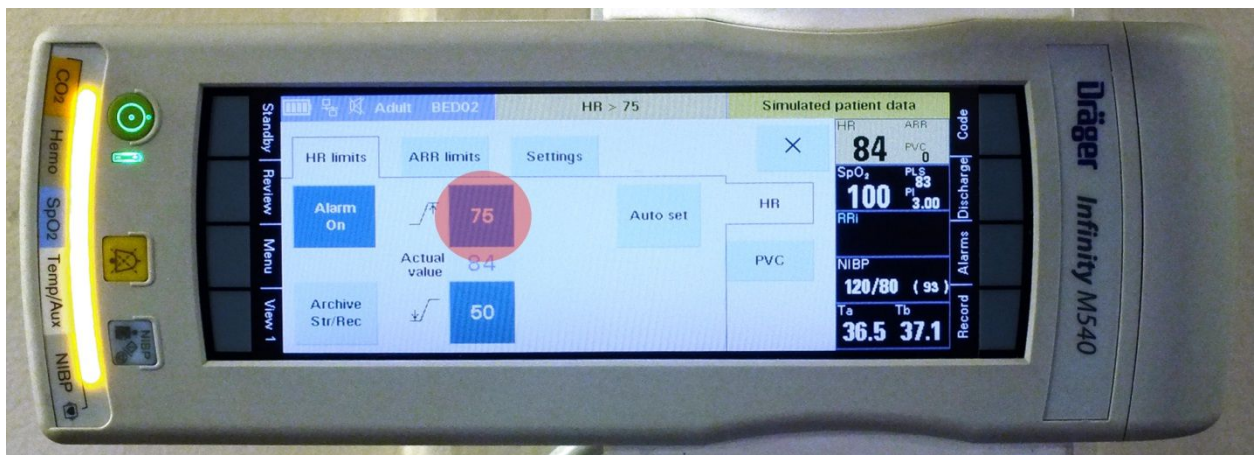


Dräger Infinity M540 Operation Manual for Demo Mode


Step 13 – Press “up arrow” button 3 times, raising upper range temp to 90.

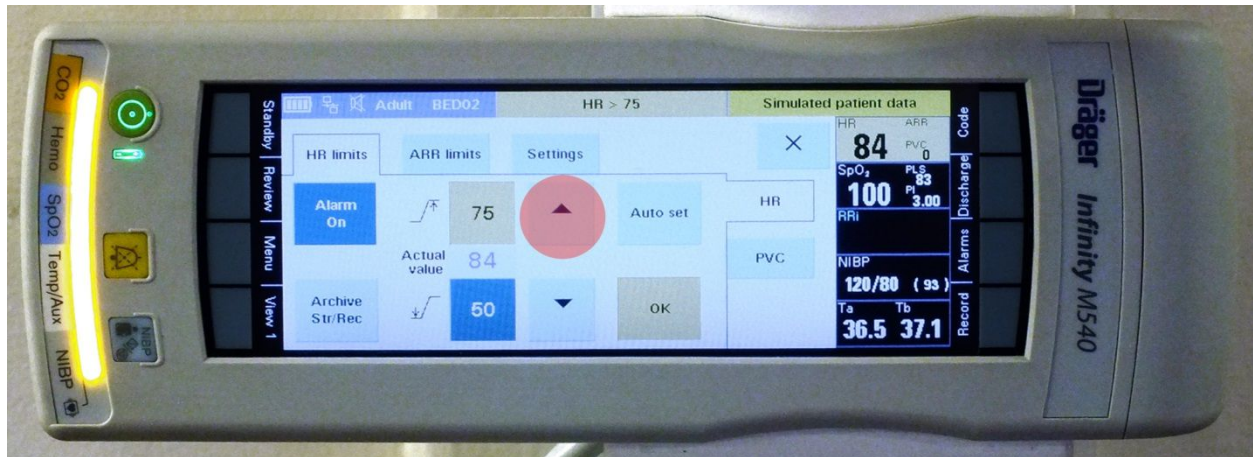

Step 14 – Press “OK” button.

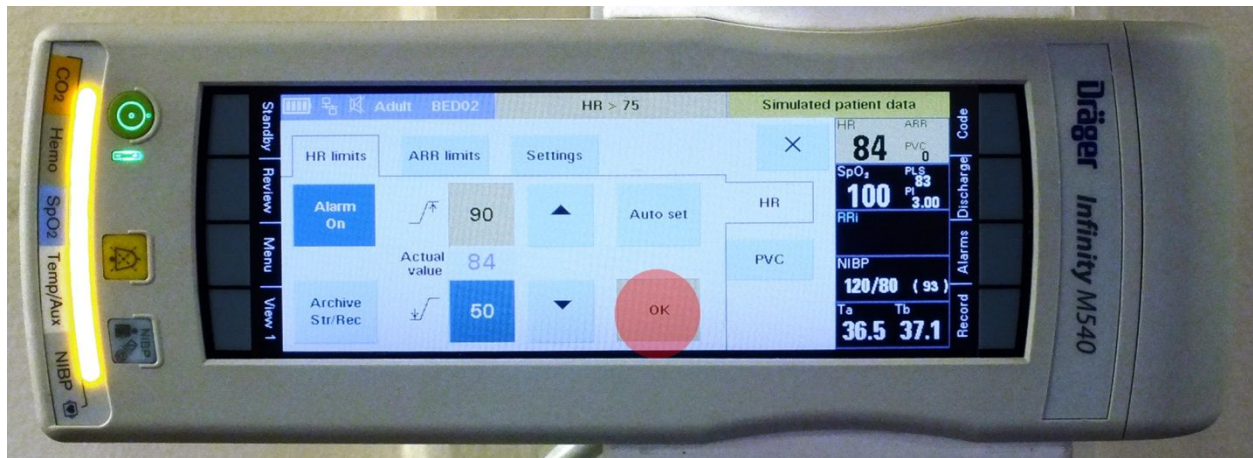


Step 15 – Press “Settings” tab.

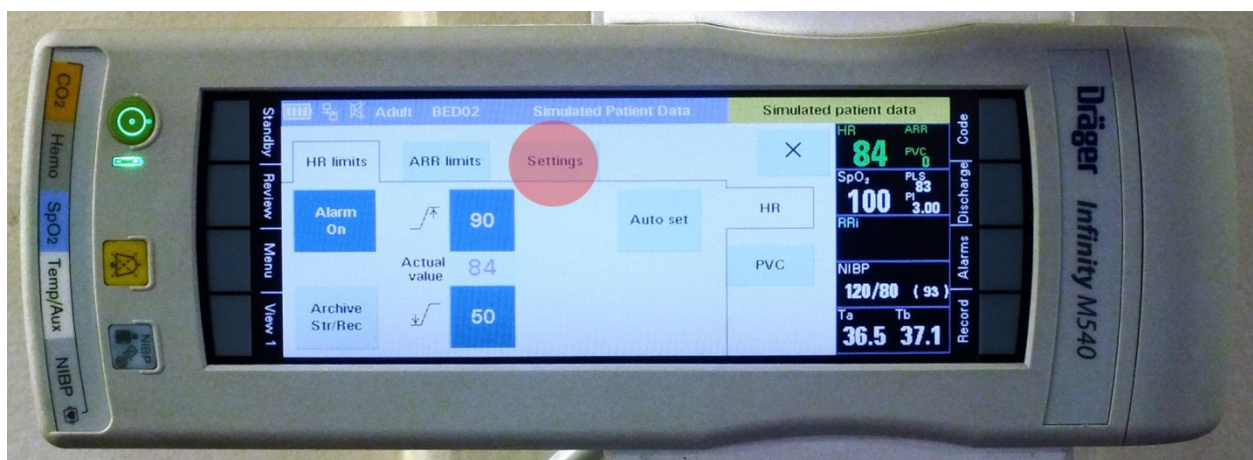


Dräger Infinity M540 Operation Manual for Demo Mode


Step 16 – Press “Tone Volume 10%” button.

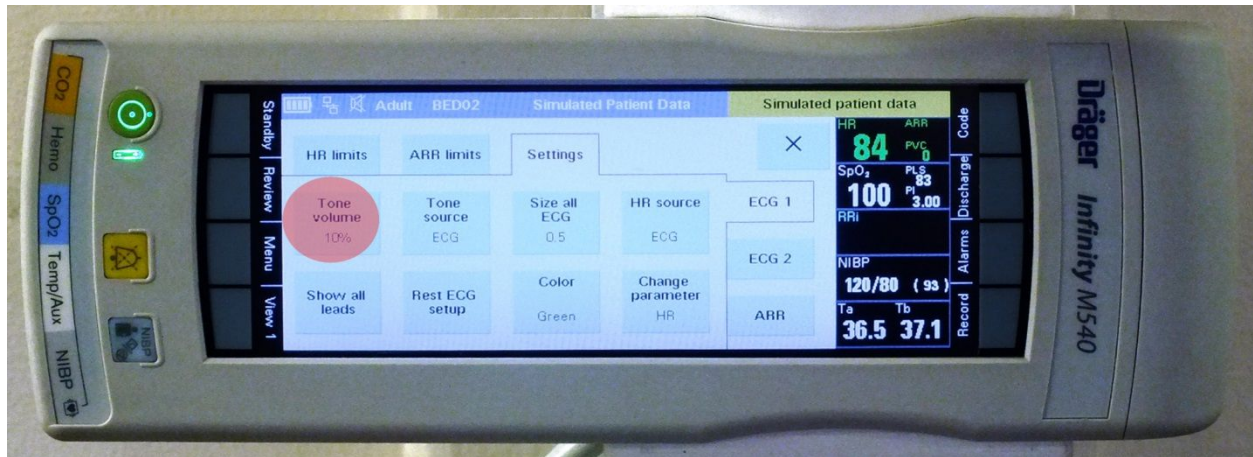

Step 17 – Press “Off” button.

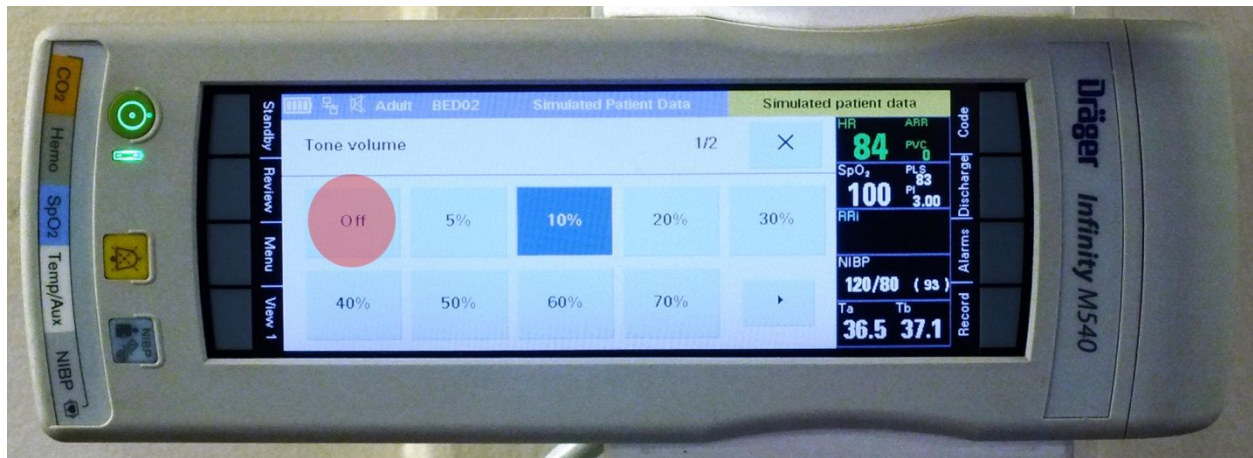

Step 18 – Press “X” to close settings screen.

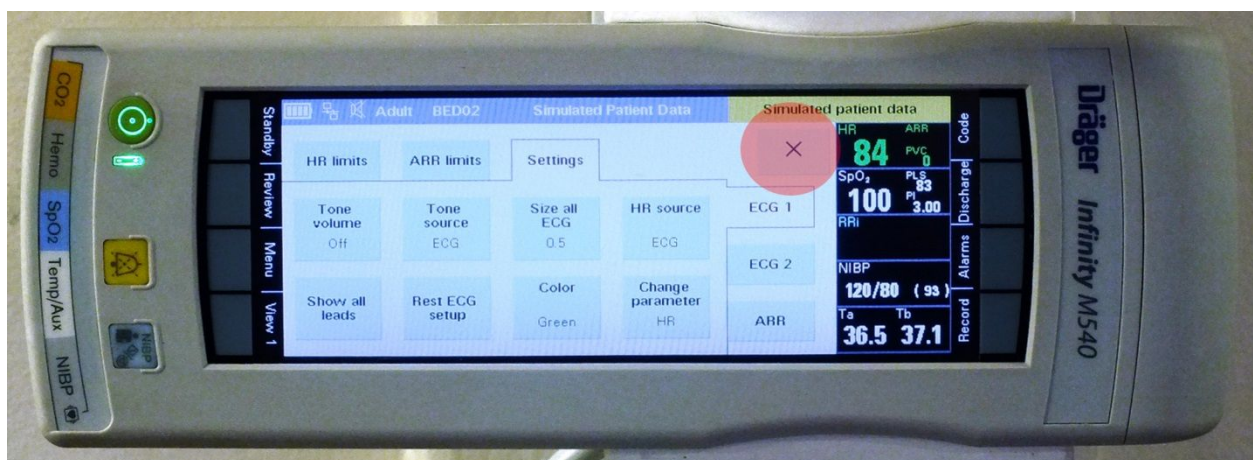


Dräger Infinity M540 Operation Manual for Demo Mode

To Shut down Dräger Infinity M540 monitor

Step 1 – Press and Hold Green Button until shutdown begins.




Step 2 – Press “Shutdown” button on touchscreen.

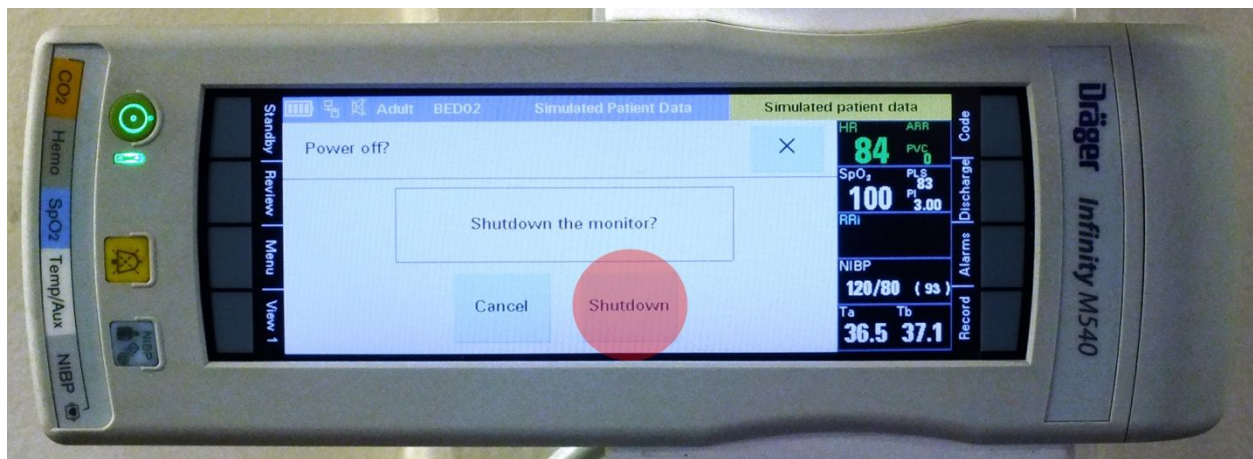

To Shut down Dräger Infinity C700 monitor

Step 1 – Press and Hold Green Button until shutdown begins.

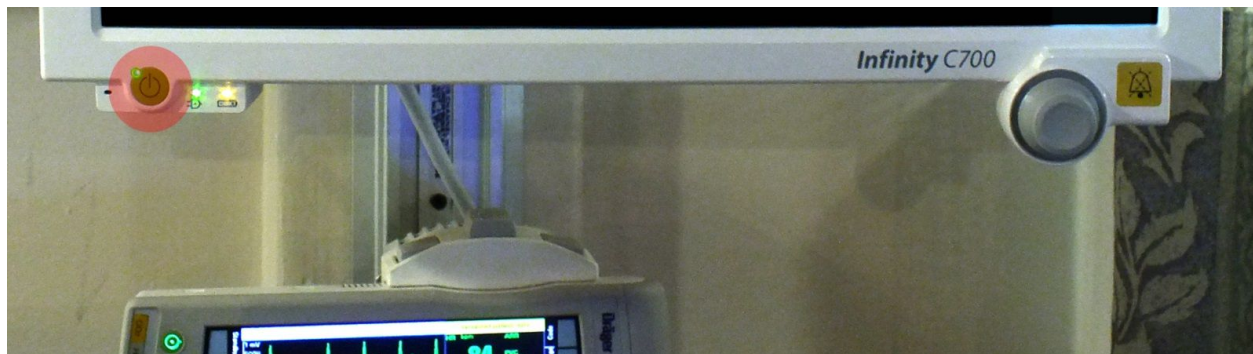

Drager Infinity M540 Operation Manual for Demo Mode

Step 2 – Press on or around “Shutdown” button on Infinity C700 monitor touchscreen.

